

GECCO ULTRAFAST FEMTOSECOND LASER DATA SHEET

SELF-STARTING FEMTOSECOND LASER

Novanta develops high-precision lasers for cutting-edge ultrafast scientific and academic research through our globally recognized brands like Laser Quantum. Designed and built for accurate results, reliable operation, and straight-forward integration, our ultrafast lasers can be found in research and development facilities around the globe.

DESIGNED FOR INTEGRATION

The gecco is a fully equipped, compact femtosecond laser in a sealed enclosure, offering an average power of >600 mW and with a <20 fs pulse duration. The gecco is offered with repetition rates of 80 MHz (70-110 MHz) with optional repetition rate locking to an external source.


Designed for simple installation and system integration, the gecco oscillator features a highly compact, vibration resistant laser head with integrated pump source, and a separate, fully featured control unit.

This highly stable laser boasts an industry leading lifetime, very low cost of ownership and comes with 2 year/5000 hours warranty covering all operating specifications.




OPTIONAL REPETITION RATE AND ACTIVE FEEDBACK

If required, active repetition rate control can be implemented for locking to an external signal source. Coarse repetition rate control is enabled by motorized mirror movement, whereas high-speed and large-amplitude feedback is realized using two piezomechanic transducers. Regulation through the TL-1000 repetition rate stabilization unit achieves a timing jitter <300 fs (0.1 Hz to 100 kHz). Alternatively, the two piezos can be driven by customer supplied electronics.

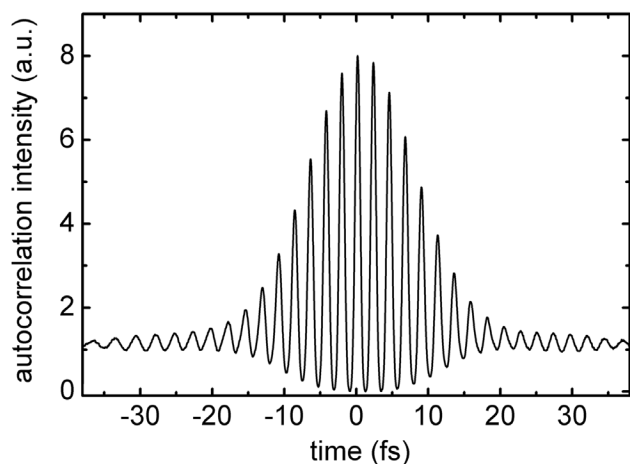

Autocorrelation trace of gecco indicating emission of pulses with 15 fs duration.

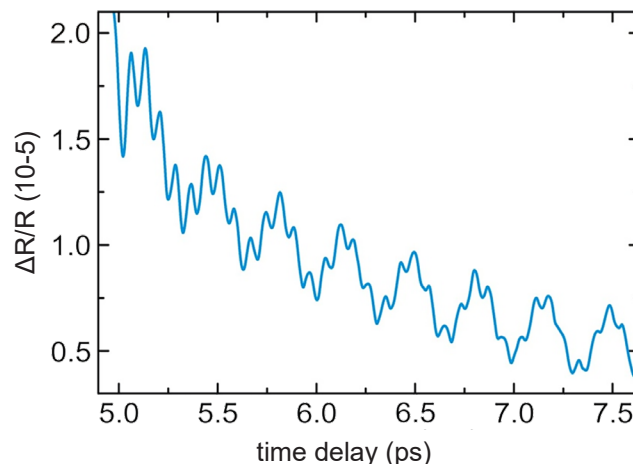

Application example: Pump-probe signature of coherent optical phonons in ZnO measured with the gecco.

GECCO ULTRAFAST FEMTOSECOND LASER SPECIFICATIONS

Specification*	gecco one
Average Power Output	> 600 mW
Center Wavelength	800 nm $\pm$ 20 nm
Pulse Duration ¹	< 20 fs
Spectral FWHM	> 40 nm
Repetition Rate ²	80 MHz (70 - 110 MHz options)
Pulse Energy ⁵	> 7.5 nJ
Beam Size ^{3,5}	0.8 mm $\pm$ 0.3 mm
Divergence ⁵	< 3 mrad
M ²	< 1.2
Power Stability (RMS within 24 hours)	< 1%
RMS Noise ⁴	< 0.1%, < 0.05% (with "Pure" option)
Polarization Ratio	> 100:1
Polarization Direction	Horizontal
Operating Temperature	21° C $\pm$ 5° C
Warm-Up Time	< 20 minutes
Weight (Head Only)	30 kg

Notes:

* Novanta operates a continuous improvement programme which can result in specifications being improved without notice.

1 - After appropriate extra cavity dispersion compensation (not included).

2 - Choose repetition rate upon order. Repetition rate accuracy $\pm$ 100 kHz.

3 - FWHM beam diameter at laser exit.


4 - Noise bandwidth 1 Hz to 1 MHz.

5 - At 80 MHz repetition rate, different at other repetition rates.

ADDITIONAL INFORMATION


- Umbilical length: 2 m
- Head weight: 15 kgs
- Cooling system included
- Warm-up time: 10 minutes
- Self-locking and maintaining
- Stable and robust
- Integrated pump laser
- The gecco is compatible with the Laser Quantum's software that allows connection to the Laser Quantum support team for monitoring laser performance and diagnosing opportunities for carrying out laser optimization.
- Pump power modulation - Modulation access to the pump power with a bandwidth of >100 kHz and modulation depth up to $\pm$ 1% is provided for feedback purposes.

DIMENSIONS (mm)

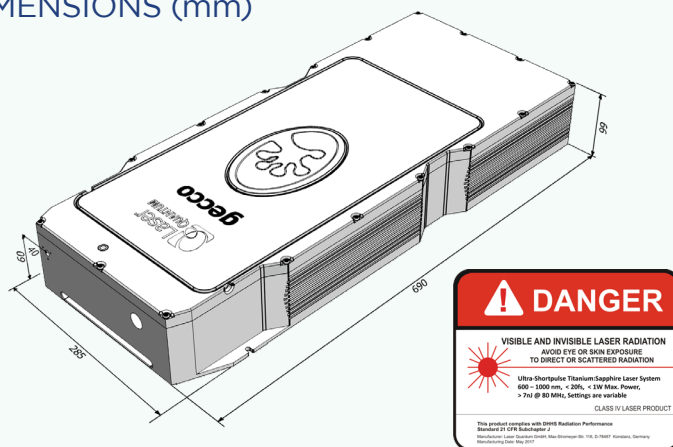

Drawings are for illustrative purposes only, please contact us for complete engineer's drawings

POWER SUPPLY UNIT

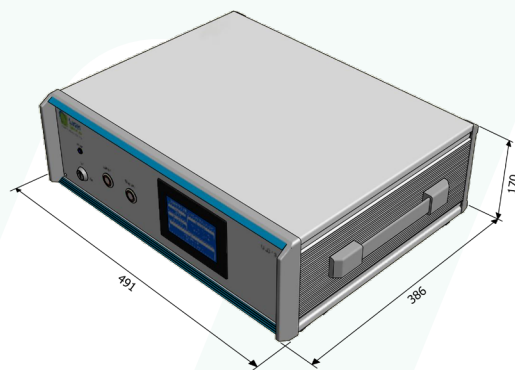

CONTACT US

Americas, Asia Pacific

Novanta Headquarters
Bedford, USA
P +1-781-266-5700

Photonics@Novanta.com

Europe, Middle East, Africa

Stockport, United Kingdom
P +44-161-975-5300

Novanta Europe GmbH
Wackersdorf, Germany
P +49 9431 7984-0

Milan, Italy
P +39-039-793-710

Photonics@Novanta.com

China

Novanta Sales & Service Office
Shenzhen, China
P +86-755-8280-5395

Suzhou, China
P +86-512-6283-7080

Photonics.China@Novanta.com


Japan

Novanta Service & Sales Office
Tokyo, Japan
P +81-3-5753-2460

Photonics.Japan@Novanta.com

光と人をつなぐ

Rayture Systems

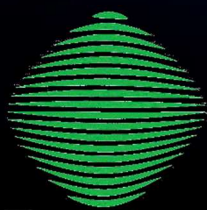

レイチャーシステムズ株式会社

〒160-0006 東京都新宿区舟町7 ロクサンビル7F

TEL : 03-3351-0717 FAX : 03-3351-6771

URL : <http://www.rayture-sys.co.jp>

E-mail : laser@rayture-sys.co.jp